## Trouble-Shooting:

Problem	Solutions
<ul style="list-style-type: none"><li>No Heat</li><li>No Power</li></ul>	<ul style="list-style-type: none"><li>Unplug unit from wall outlet.</li><li>Verify no physical damage to the cord, controller or pad.</li><li>Verify AC outlet is functional.</li><li>Restart unit by sliding control to off position then back to desired heat level.</li><li>Make sure pad is not folded or bunched in any area</li></ul>
<ul style="list-style-type: none"><li>Control Shuts Off</li></ul>	<ul style="list-style-type: none"><li>Auto off feature is set to turn unit off after two hours of use</li><li>Restart unit by sliding control to off position then back to desired heat level</li></ul>
<ul style="list-style-type: none"><li>Pad doesn't get hot enough</li></ul>	<ul style="list-style-type: none"><li>Make sure pad is not folded or bunched up in any area</li><li>Restart unit by sliding control to off position then back to desired heat level</li></ul>

**Note:** this heating pad has a built in safety feature that prevents the heating pad from being turned on immediately after the auto off feature has turned it off. Please wait at least 1 minute before turning your heating pad back on.

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## 5-YEAR (7-YEAR FOR MODEL 836) LIMITED WARRANTY

Sunbeam Products, Inc. doing business as Jarden Consumer Solutions or if in Canada, Sunbeam Corporation (Canada) Limited doing business as Jarden Consumer Solutions (collectively "JCS") warrants that for a period of five years from the date of purchase, this product will be free from defects in material and workmanship. JCS, at its option, will repair or replace this product or any component of the product found to be defective during the warranty period. Replacement will be made with a new or remanufactured product or component. If the product is no longer available, replacement may be made with a similar product of equal or greater value. This is your exclusive warranty. Do NOT attempt to repair or adjust any electrical or mechanical functions on this product. Doing so will void this warranty.

This warranty is valid for the original retail purchaser from the date of initial retail purchase and is not transferable. Keep the original sales receipt. Proof of purchase is required to obtain warranty performance. JCS dealers, service centers, or retail stores selling JCS products do not have the right to alter, modify or any way change the terms and conditions of this warranty.

This warranty does not cover normal wear of parts or damage resulting from any of the following: negligent use or misuse of the product, use on improper voltage or current, use contrary to the operating instructions, disassembly, repair or alteration by anyone other than JCS or an authorized JCS service center. Further, the warranty does not cover: Acts of God, such as fire, flood, hurricanes and tornadoes.

### What are the limits on JCS's Liability?

JCS shall not be liable for any incidental or consequential damages caused by the breach of any express, implied or statutory warranty or condition.

Except to the extent prohibited by applicable law, any implied warranty or condition of merchantability or fitness for a particular purpose is limited in duration to the duration of the above warranty.

JCS disclaims all other warranties, conditions or representations, express, implied, statutory or otherwise.

JCS shall not be liable for any damages of any kind resulting from the purchase, use or misuse of, or inability to use the product including incidental, special, consequential or similar damages or loss of profits, or for any breach of contract, fundamental or otherwise, or for any claim brought against purchaser by any other party.

Some provinces, states or jurisdictions do not allow the exclusion or limitation of incidental or consequential damages or limitations on how long an implied warranty lasts, so the above limitations or exclusion may not apply to you.

This warranty gives you specific legal rights, and you may also have other rights that vary from province to province, state to state or jurisdiction to jurisdiction

### How to Obtain Warranty Service

#### In the U.S.A. –

If you have any question regarding this warranty or would like to obtain warranty service, please call 1-800-458-8407 and a convenient service center address will be provided to you.

#### In Canada –

If you have any question regarding this warranty or would like to obtain warranty service, please call 1-800-667-8623 and a convenient service center address will be provided to you.

In the U.S.A., this warranty is offered by Sunbeam Products, Inc. doing business as Jarden Consumer Solutions located in Boca Raton, Florida 33431. In Canada, this warranty is offered by Sunbeam Corporation (Canada) Limited doing business as Jarden Consumer Solutions, located at 20 B Hereford Street, Brampton, Ontario L6Y 0M1. If you have any other problem or claim in connection with this product, please write our Consumer Service Department.

### PLEASE DO NOT RETURN THIS PRODUCT TO ANY OF THESE ADDRESSES OR TO THE PLACE OF PURCHASE.

©2008 Sunbeam Products, Inc. doing business as Jarden Consumer Solutions. All rights reserved.  
Distributed by Sunbeam Products, Inc. doing business as Jarden Consumer Solutions, Boca Raton, Florida 33431.  
1-800-458-8407. Printed in China.

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## Coussin Chauffant

## IMPORTANT NOTES

Please read this entire manual carefully before using this device.  
Product may vary from picture.

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## IMPORTANT SAFETY INSTRUCTIONS DO NOT DESTROY

**DANGER: TO REDUCE RISKS OF BURNS, ELECTRIC SHOCK, FIRE, AND PERSONAL INJURY, THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE FOLLOWING INSTRUCTIONS:**

1. READ ALL INSTRUCTIONS.
2. DO NOT USE WHILE SLEEPING.
3. DO NOT USE ON AN INFANT OR ON AN ANIMAL.
4. THIS PAD IS NOT TO BE USED ON OR BY AN INVALID, A SLEEPING OR UNCONSCIOUS PERSON, A PERSON WITH DIABETES, OR A PERSON WITH POOR BLOOD CIRCULATION.
5. DO NOT USE ON AREAS OF INSENSITIVE SKIN.
6. BURNS CAN OCCUR REGARDLESS OF CONTROL SETTING. CHECK SKIN UNDER PAD FREQUENTLY.
7. DO NOT USE IN AN OXYGEN ATMOSPHERE.
8. DO NOT USE PINS OR OTHER METALLIC MEANS TO FASTEN THIS PAD IN PLACE.
9. DO NOT SIT ON OR AGAINST OR CRUSH PAD-AVOID SHARP FOLDS.
10. NEVER PULL THIS PAD BY THE SUPPLY CORD.
11. DO NOT USE THE CORD AS A HANDLE.
12. CAREFULLY EXAMINE INNER COVER BEFORE EACH USE. DISCARD THE PAD IF INNER COVERING SHOWS ANY SIGNS OF DETERIORATION.
13. NEVER USE PAD WITHOUT THE REMOVABLE COVER IN PLACE.
14. THIS UNIT SHOULD NOT BE USED BY CHILDREN WITHOUT ADULT SUPERVISION.
15. PLACE PAD ON TOP OF AND NOT UNDER THE PART OF THE BODY NEEDING HEAT.
16. USE THIS PAD ONLY ON A 110-120 VOLT AC CIRCUIT. UNPLUG WHEN NOT IN USE.
17. DO NOT TAMPER WITH THIS PAD IN ANY WAY.
18. THERE ARE NO USER-SERVICEABLE PARTS. IF FOR ANY REASON THIS PAD DOES NOT FUNCTION SATISFACTORILY, SEE WARRANTY FOR CONSUMER SERVICE INFORMATION.
19. DO NOT USE THIS PAD WITH LINIMENT, SALVE OR OINTMENT PREPARATIONS, ESPECIALLY ONES THAT CONTAIN HEAT- PRODUCING INGREDIENTS. SKIN BURNS COULD RESULT.
20. DO NOT BEND OR PINCH CORD AND DO NOT USE IF CORD IS DAMAGED.
21. NEVER LEAVE THE APPLIANCE UNATTENDED, ESPECIALLY IF CHILDREN ARE PRESENT.
22. LOOP CORD LOOSELY WHEN STORING. TIGHT WRAPPING MAY DAMAGE CORD AND INTERNAL PARTS.
23. SAVE THESE INSTRUCTIONS.

This appliance has a polarized plug (one blade is wider than the other). This plug fits a polarized outlet only one way; if the plug does not fit fully into the outlet, reverse the plug. If it still does not fit, contact a qualified electrician.

DO NOT ATTEMPT TO DEFEAT THIS SAFETY FEATURE.

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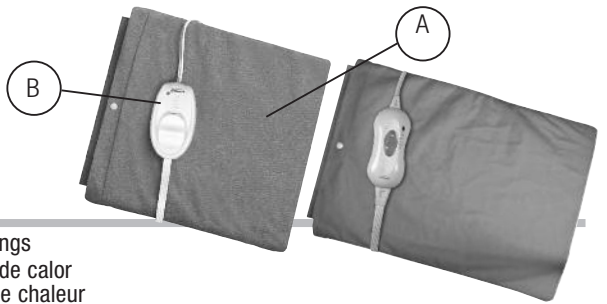
### HOW TO USE YOUR HEATING PAD

1. Lay the cloth cover on a flat surface with the printing face up. Open the snap at the bottom.
2. Insert the pad into the cover completely, fitting one side at a time. (Figure 1) Make sure the cord and control remain on the outside of the cover.
3. The pad will lie flat, filling the entire cover. Close snap.
4. Plug into any 110-120 volt AC electrical outlet and set control to the desired heat.

**WARNING:** Never use the heating pad without the cloth cover.

### FOR MODELS WITH HOT (RED) GEL PACK:

1. Remove gel pack from heating pad.
  - a. Place gel pack in microwave.
  - b. Make sure gel is spread evenly across pack for proper heating

### CAUTION: UNEVEN GEL COULD CAUSE HOT SPOTS

- c. Heat at full power (700 Watts) for one minute. NOTE: Heat time may vary based on microwave wattage.
- d. Remove and check heat guide.
2. If not hot enough, return to microwave and heat in ten second intervals, checking Heat Guide until ready.
3. If Heat Guide reads TOO HOT, let cool before use.

### GEL PACK SAFETY REMINDERS:

- a. Gel Pack can be used without heating pad, but must be enclosed in cover.
- b. NEVER HEAT GEL PACK FOR MORE THAN 90 SECONDS AT ONE TIME OR GEL PACK MAY RUPTURE AND LEAK.
- c. NEVER BOIL GEL PACK.

4. Reassemble gel pack to heating pad by aligning hook and loop with the heating pad.

### Insert Heating Pad and Gel Pack in cover.

- a. Lay heating pad flat.
- b. Attach Gel Pack to heating pad. NOTE: You can use heating pad WITHOUT Gel pack.
- c. Slip inside cover, making sure cord and controls remain on outside.
- d. Close snap.

6. Plug cord in wall. NOTE: Although your heating pad is meant to get hot, it should NEVER hurt or feel like it is burning your skin. If it does, remove the heating pad immediately.

### FOR MODELS WITH COLD (BLUE) GEL PACK:

Follow the instructions located on the gel pack.

### FOR MODELS WITH AUTO OFF FEATURES

The Auto-off feature is indicated by a LED light in the upper left corner of the controller. The LED light is red when the heating pad is turned on. After 2 hours of continuous usage, the pad shuts off automatically and the LED light turns off. **Note: this heating pad has a built in safety feature that prevents the heating pad from being turned on immediately after the auto off feature has turned it off. Please wait at least 1 minute before turning your heating pad back on.**

### USING MOIST HEAT OPTION ON SELECTED MODELS

Your heating pad may include a separate sponge for moist heat.

1. Run warm water over sponge until it is thoroughly saturated.
2. Gently wring out excess water.
3. Place the sponge on top of the heating pad and slip the cloth cover over the heating pad. Snap cover closed.

**DANGER:** Do not immerse pad in water or place under running water to dampen surface. Failure to comply may result in personal injury.

### FABRIC COVER WASHING INSTRUCTIONS

1. Disconnect pad from outlet.
2. Remove cover from pad. (Figure 1)
3. Machine or hand wash; use mild detergent. Wash with cool water only.
4. Hang to dry or machine dry on gentle, cool cycle. (Figure 2)

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