Renue® Heat Wrap with XpressHeat®

Model: 885

IMPORTANT SAFETY INSTRUCTIONS DO NOT DESTROY

DANGER: TO REDUCE THE RISK OF BURNS, ELECTRIC SHOCK AND FIRE, THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE FOLLOWING INSTRUCTIONS:

1. READ ALL INSTRUCTIONS.
2. DO NOT USE WHILE SLEEPING.
3. DO NOT USE ON AN INFANT OR ON AN ANIMAL.
4. THIS PAD IS NOT TO BE USED ON OR BY AN INVALID, SLEEPING, HELPLESS OR UNCONSCIOUS PERSON, A PERSON WITH DIABETES, A PERSON WITH POOR BLOOD CIRCULATION OR INSENSITIVITY TO HEAT, OR A PERSON INCAPACITATED BY MEDICATION.
5. THIS APPLIANCE IS NOT TO BE USED BY CHILDREN. CLOSE SUPERVISION IS NECESSARY WHEN ANY APPLIANCE IS USED NEAR CHILDREN. CHILDREN SHOULD NOT PLAY WITH THE APPLIANCE.
6. DO NOT USE ON AREAS OF INSENSITIVE SKIN.
7. BURNS CAN OCCUR REGARDLESS OF CONTROL SETTING, CHECK SKIN UNDER PAD FREQUENTLY.
8. IF PROVIDED WITH A REMOVEABLE COVER, NEVER USE PAD WITHOUT COVER IN PLACE.
9. DO NOT USE IN AN OXYGEN ATMOSPHERE.
10. DO NOT USE PINS OR OTHER METALLIC MEANS TO FASTEN THIS PAD IN PLACE.
11. DO NOT SIT ON, AGAINST, OR CRUSH PAD – AVOID SHARP FOLDS. PLACE PAD ON TOP OF AND NOT UNDER THE PART OF BODY NEEDING HEAT.
12. UNPLUG WHEN NOT IN USE. NEVER LEAVE APPLIANCE UNATTENDED, ESPECIALLY IF CHILDREN ARE PRESENT.
13. NEVER PULL THIS PAD BY THE SUPPLY CORD. DO NOT USE THE CORD AS A HANDLE. DO NOT PINCH CORD.


14. DO NOT TUCK IN, TRAP, FOLD, CROSS, PINCH, OR WRAP CONTROL OR CORD WITHIN HEATING PAD OR ANY OTHER INSULATING MATERIAL DURING USE.
15. DO NOT PLACE ANYTHING OVER THE HEATING PAD WHILE IN USE OR PLUGGED IN, SUCH AS A TOWEL, BLANKET, CLOTHING, OR ANY OTHER INSULATING MATERIAL, AND NEVER ALLOW HEATING PAD TO WRAP AROUND ITSELF.
16. CAREFULLY EXAMINE PAD, COVER, CORD, AND CONTROLLER BEFORE EACH USE. DISCARD THE PAD IF ANY SIGN OF DETERIORATION, WEAR OR DAMAGE. DO NOT TAMPER WITH THIS PAD IN ANY WAY. THERE ARE NO USER SERVICABLE PARTS. IF PAD DOES NOT FUNCTION SATISFACTORILY, SEE WARRANTY FOR CUSTOMER CONTACT INFORMATION.
17. LOOP CORD LOOSELY WHEN STORING, TIGHT WRAPPING MAY DAMAGE CORD AND INTERNAL PARTS. ALLOW TO COOL BEFORE STORING AND DO NOT CREASE BY PLACING ITEMS ON TOP OF IT. STORE IN THE ORIGINAL PACKAGING, IN A COOL, DARK, DRY PLACE.
18. IF PROVIDED WITH A DETACHABLE CONTROLLER, USE ONLY WITH THE PROVIDED CONTROLLER.
19. USE THIS PAD ONLY AT RATED VOLTAGE, AS MARKED ON THE PRODUCT. DO NOT USE WITH GENERATORS, POWER CONVERTERS, OR INVERTERS.
20. DO NOT USE THIS PAD WITH LINIMENT, SALVE OR OINTMENT PREPERATIONS. SKIN BURNS COULD RESULT.
21. THIS APPLIANCE IS NOT INTENDED FOR MEDICAL USE IN HOSPITALS.
22. SOME MODELS MAY CONTAIN FILLING MATERIAL (PELLETS OR BEADS), MADE FROM GLASS, ABS PLASTIC (NON-MICROWAVEABLE MODELS) OR SILICON DIOXIDE (MICROWAVEABLE MODELS). IF FILLING MATERIAL IS RELEASED DUE TO PRODUCT DAMAGE, DISCONTINUE USE OF PRODUCT. DO NOT INGEST FILLING MATERIAL. IF FILLING MATERIAL IS INGESTED, SEEK MEDICAL ATTENTION.
23. SAVE THESE INSTRUCTIONS.

HOUSEHOLD USE ONLY

NORTH AMERICAN MODELS WITH POLARIZED PLUGS:


POLARIZED PLUG

THIS APPLIANCE HAS A POLARIZED PLUG (ONE BLADE IS WIDER THAN THE OTHER). AS A SAFETY FEATURE TO REDUCE THE RISK OF ELECTRICAL SHOCK, THIS PLUG IS INTENDED TO FIT IN A POLARIZED OUTLET ONLY ONE WAY. IF THE PLUG DOES NOT FIT FULLY IN THE OUTLET, REVERSE THE PLUG. IF IT STILL DOES NOT FIT, CONTACT A QUALIFIED ELECTRICIAN. DO NOT ATTEMPT TO DEFEAT THIS SAFETY FEATURE OR MODIFY THE PLUG IN ANY WAY. IF THE PLUG FITS LOOSELY INTO THE AC OUTLET OR IF THE AC OUTLET FEELS WARM DO NOT USE THAT OUTLET.

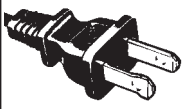Renue® Heat Wrap is specifically designed to fit the neck and shoulder area, where users experience many physical signs of tension and stress. Drape the Renue® Heat Wrap comfortably around your shoulders making sure that the neck flap is properly in place. Weighted edges and a hook closure help to secure Renue® Heat Wrap in place.

Figure 1

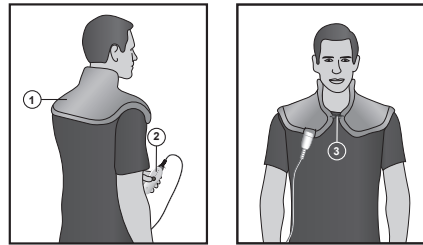HEATING PAD COMPONENTS

1. Exterior Pad Fabric
2. Controller
3. Closure

HOW TO APPLY YOUR HEATING PAD

Heating pads are effective pain relievers to soothe muscle and joint pain. The Renue® Heat Wrap with XpressHeat® Technology is specifically designed to fit the back of the neck, shoulders, and upper back.

- NOTE: Always consult with your physician or health professional for recommended treatment of an injury or chronic muscle or joint pain.
1. Drape the pad comfortably over your shoulders.
 2. Once in place, secure the pad by connecting the hook closures on either end of the shoulders (Figure 1, Number 3).
 3. To remove pad, detach the hook closures, and lift the pad off of your shoulders.

HOW TO USE YOUR HEATING PAD

1. Make sure the Control Connector (Figure 3, Letter A) is securely attached to the Pad Connector (Figure 3, Letter B) until they are tightly joined together.
2. Plug into 120V AC electrical outlet.
3. To turn the pad on, press the Heat Control Button (Figure 2, Letter A).
4. Select from 4 heat settings (High, Medium, Low, and Warm) by repeatedly pressing the Heat Control Button until the desired heat setting indicator lights up (Figure 2, Letter B).
5. To turn the pad off, press the Heat Control Button repeatedly until all of the indicator lights turn off.

CAUTION:

- DO NOT SIT ON, LEAN AGAINST, CRUSH OR FOLD PAD DURING USE.
- DO NOT USE WHILE SLEEPING.
- SEE ADDITIONAL WARNINGS IN THE IMPORTANT SAFETY INSTRUCTIONS SECTION.

AUTO-OFF FEATURE

Your heating pad controller includes an Auto-Off feature. This feature is designed to conserve energy by automatically turning off after 2 hours of continuous use. The Auto-Off Indicator (Figure 2, Letter C) will flash when the Auto-Off feature is activated and the pad shuts off. To reset after Auto-Off, press the Heat Control Button until the desired heat setting is reached.

NOTE: After a power outage, the Auto-Off feature automatically resets the time.

HOW TO USE YOUR HEATING PAD WITH MOIST HEAT THERAPY

Moist heat provides deep, penetrating pain relief. When using moist heat therapy, you should protect furniture and clothing from becoming damp.

1. Dampen the body side of the heating pad with water by patting the surface with a wet sponge or cloth, or use a spray bottle to apply a fine mist. Do not thoroughly soak pad.
2. Do not spray or dampen the pad connector or the control connector.

DANGER: DO NOT IMMERSER THE PAD OR WIRING IN WATER OR PLACE UNDER RUNNING WATER. FAILURE TO COMPLY MAY RESULT IN PRODUCT FAILURE AND/OR PERSONAL INJURY. THE ELECTRICAL PAD CONNECTOR MUST BE THOROUGHLY DRY BEFORE ESTABLISHING ELECTRICAL CONNECTION.

WASHING INSTRUCTIONS

CAUTION: DO NOT DRY CLEAN! DO NOT BLEACH! DO NOT USE WRINGER! DO NOT IRON!

1. Disconnect pad from outlet. DO NOT PULL BY THE CORD.
2. Disconnect controller from pad by firmly pulling Control Connector (Figure 3, Letter A) and the Pad Connector (Figure 3, Letter B) away from each other until they are completely separated.
3. Machine wash pad in cold water on gentle cycle or hand wash. Tumble dry on low heat and remove promptly from dryer or drape to dry.
4. DO NOT plug Control Connector into Pad Connector until completely dry.
5. For information on service, visit our website at www.Sunbeam.com.

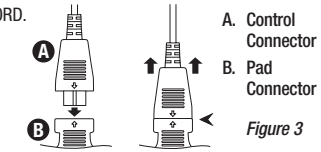
Figure 3

TROUBLE SHOOTING

Problem	Solution
<ul style="list-style-type: none"> • No Heat • No Power 	<ol style="list-style-type: none"> 1. Verify tight connection between pad and controller. 2. Verify AC outlet is functional. 3. Verify no physical damage to the cord, controller or pad. If failure continues, discontinue use. See warranty information to contact Sunbeam Consumer Service.
<ul style="list-style-type: none"> • Pad Not Hot Enough 	<ul style="list-style-type: none"> • Ensure that the body side of the pad is laying against the body.
<ul style="list-style-type: none"> • Pad has No Heat & Controller 'Auto-Off' Indicator is Flashing 	Auto-Off feature is set to turn the unit off after 2 hours of continuous use. <ul style="list-style-type: none"> • To re-start unit, press the Heat Control Button until desired heat setting is reached.
<ul style="list-style-type: none"> • Pad has a Hot Spot 	<ul style="list-style-type: none"> • Ensure the pad is not folded, bunched, or crushed in any area.
<ul style="list-style-type: none"> • 'High' Heat Setting Indicator is Flashing 	The 'High' Heat Setting Indicator light will start flashing when connection is broken between pad and controller while using the product. <ol style="list-style-type: none"> 1. To reset, simply connect pad & controller together and press the Heat Control Button until desired heat setting is reached. 2. If flashing light persists, an error has occurred. Discontinue use of the product and see warranty information to contact Sunbeam Customer Service.


FCC CERTIFIED

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

NOTE: This product has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This product generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this product does cause harmful interference to radio or television reception, which can be determined by turning the product off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
 - Increase the separation between the product and receiver.
 - Connect the product into an outlet on a circuit different from that to which the receiver is connected.
 - Consult the dealer or an experienced radio/TV technician for help.
- Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Figure 2

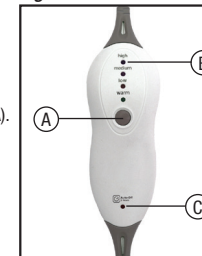- A. Heat Control Button
- B. Heat Setting Indicator
- C. Auto-Off Indicator

5 YEAR LIMITED WARRANTY

Sunbeam Products, Inc. ("Sunbeam") or if in Canada, Newell Brands Canada ULC ("NBC") warrants that for a period of five (5) years from the date of initial retail purchase, this product will be free from defects in material and workmanship, when used in accordance with the instructions provided.

Sunbeam (or NBC), at its option, will repair or replace this product or any component of the product found to be defective during the warranty period.

Replacement will be made with a new or remanufactured product or component. If the product is no longer available, replacement may be made with a similar product of equal or greater value. This is your exclusive warranty. Do NOT attempt to repair or adjust any electrical or mechanical functions on this product. Doing so will void this warranty.

This warranty is valid for the original retail purchaser from the date of initial retail purchase and is not transferable. Keep the original sales receipt. Proof of purchase is required to obtain warranty performance. Sunbeam (or NBC) dealers, service centers, or retail stores selling Sunbeam (or NBC) products do not have the right to alter, modify or any way change the terms and conditions of this warranty.

This warranty does not cover normal wear of parts or damage resulting from any of the following: negligent use or misuse of the product, use on improper voltage or current, use contrary to the operating instructions, disassembly, repair or alteration by anyone other than Sunbeam (or NBC) or an authorized Sunbeam (or NBC) service center. Further, the warranty does not cover: Acts of God, such as fire, flood, hurricanes and tornadoes.

What are the limits on Sunbeam's (or NBC's) Liability?

Sunbeam (or NBC) shall not be liable for any incidental or consequential damages caused by the breach of any express, implied or statutory warranty or condition. Except to the extent prohibited by applicable law, any implied warranty or condition of merchantability or fitness for a particular purpose is limited in duration to the duration of the above warranty. Sunbeam (or NBC) disclaims all other warranties, conditions or representations, express, implied, statutory or otherwise. Sunbeam (or NBC) shall not be liable for any damages of any kind resulting from the purchase, use or misuse of, or inability to use the product including incidental, special, consequential or similar damages or loss of profits, or for any breach of contract, fundamental or otherwise, or for any claim brought against purchaser by any other party. Some provinces, states or jurisdictions do not allow the exclusion or limitation of incidental or consequential damages or limitations on how long an implied warranty lasts, so the above limitations or exclusion may not apply to you. This warranty gives you specific legal rights, and you may also have other rights that vary from province to province, state to state or jurisdiction to jurisdiction.

How to Obtain Warranty Service

In the U.S.A.

If you have any question regarding this warranty or would like to obtain warranty service, please call 1 800 435-1250 and a convenient service center address will be provided to you.

In Canada

If you have any question regarding this warranty or would like to obtain warranty service, please call 1 800 667- 8623 and a convenient service center address will be provided to you.

In the U.S.A., this warranty is offered by Sunbeam Products, Inc. located in Boca Raton, Florida 33431.

In Canada, this warranty is offered by Newell Industries Canada ULC, located at 20B Hereford Street, Brampton, Ontario L6Y 0M1. If you have any other problem or claim in connection with this product, please write our Consumer Service Department.

PLEASE DO NOT RETURN THIS PRODUCT TO ANY OF THESE ADDRESSES OR TO THE PLACE OF PURCHASE.

For inquiries regarding recycling and proper disposal of this product, please contact your local waste management facility.

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